


Q-Swab™

Environmental Swab-Collection System

Q-Swab is a ready-to-use dilution and delivery device to be used in the food and beverage industry for environmental surface sampling. Money and time are saved by eliminating the need for pipettes and broth preparation. Each of the four broths offered in Q-Swab are self-contained in a 1.0 mL reservoir by patented Snap Valve™ technology. When the Snap Valve is snapped, broth is released down the swab shaft, neutralizing residual sanitizers and facilitating recovery of bacteria. The sample is then ready to be poured onto dehydrated media or streaked onto Agar plates for culturing. Q-Swab can be used either wet or dry to sample environmental surfaces. Simply swab, snap and squeeze, and a sample is ready for transport to the laboratory.

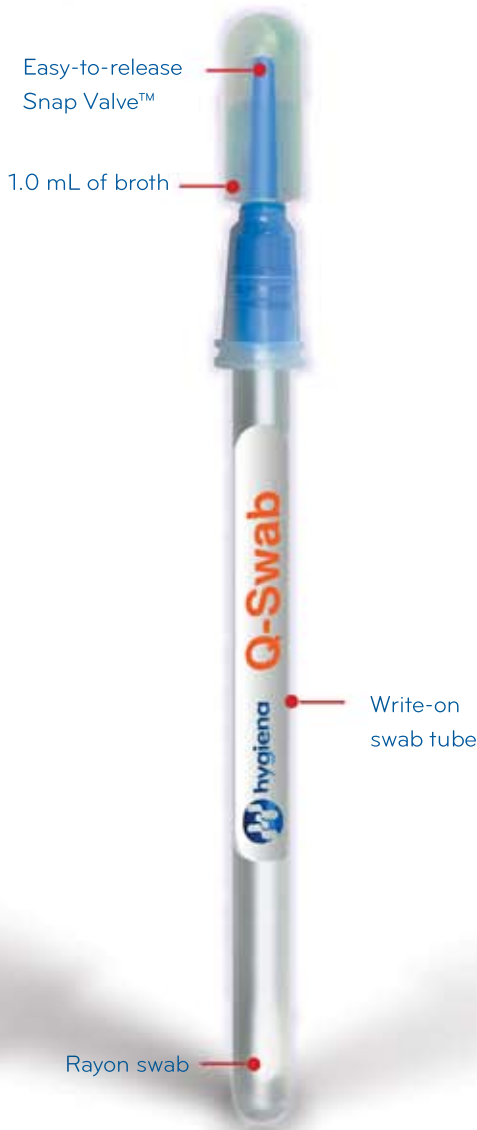

Key Features	
✔	Aseptically filled and gamma irradiated
✔	All-in-one sample collection system
✔	12-month shelf life
✔	1.0 mL fill of broth - perfectly hydrates dried media
✔	Easy to use
Ideal for	
✔	Food and beverage industry



Americas
 941 Avenida Acaso
 Camarillo, CA 93012
 USA
Tel: 805.388.8007
Fax: 805.388.5531
Email: enquiries@hygiena.net
Web: www.hygiena.net

Europe
 Unit 11 Wenta Business Centre
 Colne Way Watford
 Hertfordshire WD2 4ND UK
Tel: +44 (0)1923 818821
Fax: +44 (0)1923 818825
Email: enquiries@hygiena.net
Web: www.hygiena.net

ORDER #	
QS-1000	Buffered Peptone Water
QS-1200	Lethen Broth
QS-1300	Stuart's Medium
QS-1400	Amies Medium

(250 tests per box)